

Download the file

If you are using Internet Explorer in Windows 7, you can either right-click on the zip file, or left click on it.

If you left click on the file, you should get a pop up with a prompt to open or save the file. Please save the file, first. Then, Extract all files from the saved zip file on your local computer.

Left-click download

The screenshot shows a web page with a sidebar on the left containing links for 'Industry Resources', 'Web Seminars', 'Spectrum Viewer', and 'Software Licensing', along with a 'Contact EDAX' button. The main content area features a link for 'Spectrum Viewer Win 7 Zip' (485 KB) with a red box around it. Below the link is a description of the software. At the bottom of the page, there are social media links for LinkedIn and YouTube, and a navigation menu. A Windows 7 download dialog box is overlaid at the bottom, asking 'Do you want to open or save Spectrum_Viewer_Windows_7.zip from edax.com?'. The 'Save' button in the dialog is highlighted with a red box, and a red arrow points from the 'Spectrum Viewer Win 7 Zip' link to this button.


Right-click download

If you Right click on the file, choose Save Target As....

Install the software

Then, you can run the install program from the Extracted files. Be sure you have Administrator rights.

After you extract the files, right-click on **InstallToC.bat** and select to **Run as Administrator**. You should get three prompts showing a process succeeded.

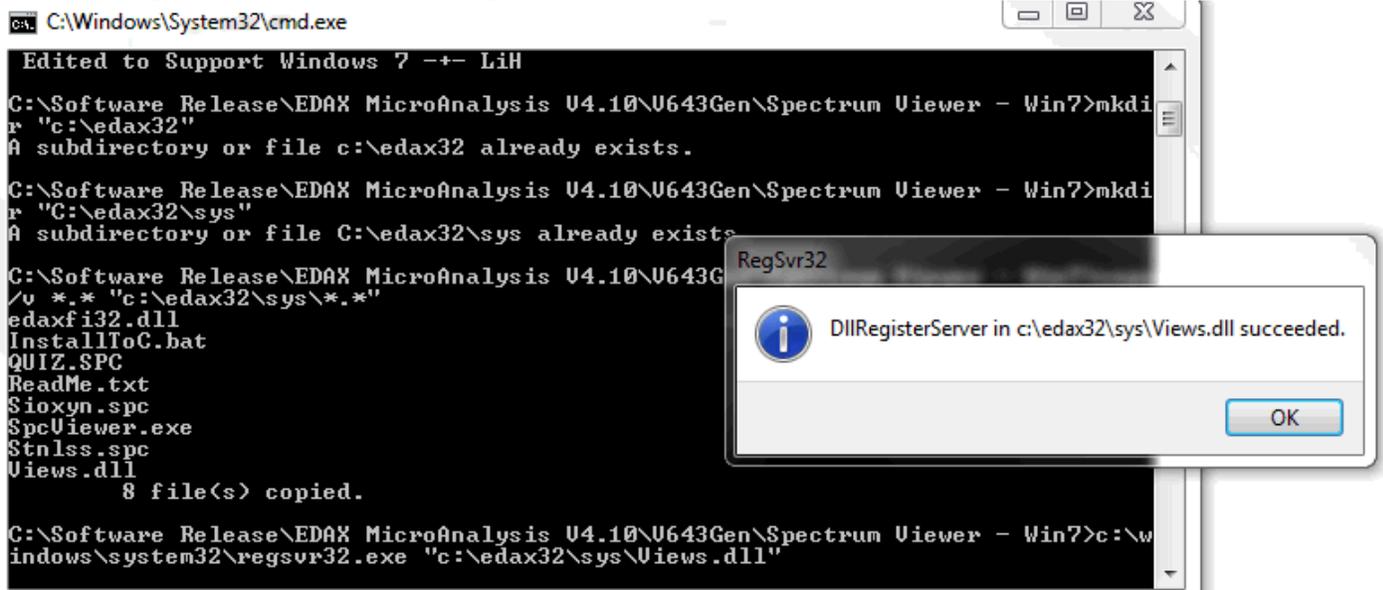

```
C:\Windows\System32\cmd.exe
Edited to Support Windows 7 -- LiH
C:\Software Release\EDAX MicroAnalysis U4.10\U643Gen\Spectrum Viewer - Win7>mkdir
r "c:\edax32"
A subdirectory or file c:\edax32 already exists.
C:\Software Release\EDAX MicroAnalysis U4.10\U643Gen\Spectrum Viewer - Win7>mkdir
r "C:\edax32\sys"
A subdirectory or file C:\edax32\sys already exists
C:\Software Release\EDAX MicroAnalysis U4.10\U643Gen\Spectrum Viewer - Win7>copy
/u *.* "c:\edax32\sys\*.*"
edaxfi32.dll
InstallToC.bat
QUIZ.SPC
ReadMe.txt
Sioxyn.spc
SpcViewer.exe
Stnlss.spc
Views.dll
      8 file(s) copied.
C:\Software Release\EDAX MicroAnalysis U4.10\U643Gen\Spectrum Viewer - Win7>c:\w
indows\system32\regsvr32.exe "c:\edax32\sys\Views.dll"
```

The program will be installed in **c:\edax32\sys** folder.
Run **SpcViewer** from that folder.